

Troubleshooting

If problem cannot be corrected, consult your Dealer/Home Care equipment provider.

PROBLEM	PROBABLE CAUSE	REMEDY
A. No flow	1. Cylinder valve closed	1. Turn on cylinder
	2. Conserver in "OFF" position	2. Set to prescribed setting
	3. Cylinder empty	3. Replace cylinder
	4. Dial not at setting	4. Set dial so indicator points to a setting
	5. Conserver not sensing breath	5. Check position of cannula in nose Do not breathe through mouth
B. Leaking at Cylinder Connection	1. Missing or defective O-ring Seal	1. Replace O-ring Seal

Precision Medical, Inc.
300 Held Drive • Northampton, PA 18067
800-272-7285 • 610-262-6090 • fax: 610-262-6080

504442 ?M 4/05 Printed in USA

Setting up your

EASYPULSE5

Oxygen Conserving Regulator

1900 Series (CGA 540 Connection)



PRECISION MEDICAL®




Quick Start Guide only

Refer to User Manual for complete operating instructions, warnings and cautions.

EasyPulse5 Oxygen Conserving Regulator (Conserver)


PARTS

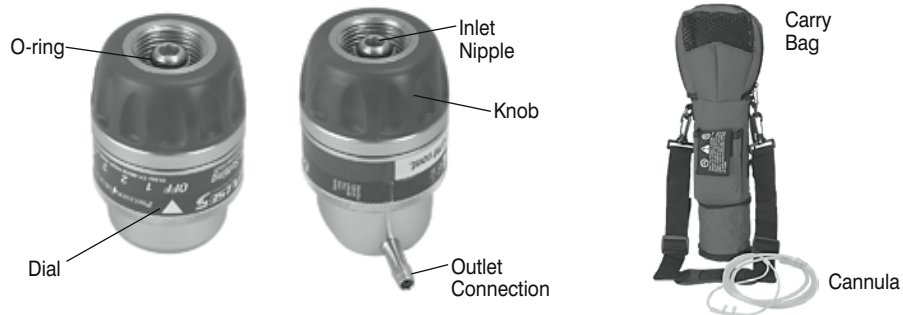

STEP 1

Inspect Conserver and Cylinder Valve for damage or contaminants. Do NOT use if any are found. “Crack” cylinder valve by opening the valve slightly and closing.

Cylinder Valve

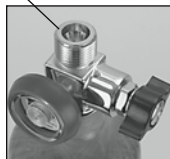

STEP 2

Ensure Conserver is in the OFF position and O-ring Seal is in place. Align Nipple with recessed hole in the cylinder. Hand tighten Knob. Do NOT use tools.

View of properly attached conserver




STEP 3

When using a Precision Medical Carry Bag, slide the cylinder into the bag with the gauge facing the mesh. Attach the cannula to the Conserver’s Outlet Connection.



STEP 4

Slowly open the cylinder valve. Turn Dial to the prescribed setting. Zip the bag closed.

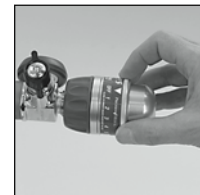

STEP 5

Place cannula in nostrils and breathe through nose to begin Oxygen flow. If using a bag adjust the straps for maximum comfort.



Cannula prongs in nostrils



Typical bag carrying configurations

STEP 6

When finished using: Turn the cylinder valve to the OFF position. Turn Conserver Dial to “2 LPM Continuous” position. Wait for Oxygen to stop flowing from Conserver. Remove Conserver from Cylinder Valve. Store Conserver and Carry Bag in a safe place.